
ICONNECT

NAVISTAR CABLE

INSTRUCTIONS FOR USE

REF 001-009070-1

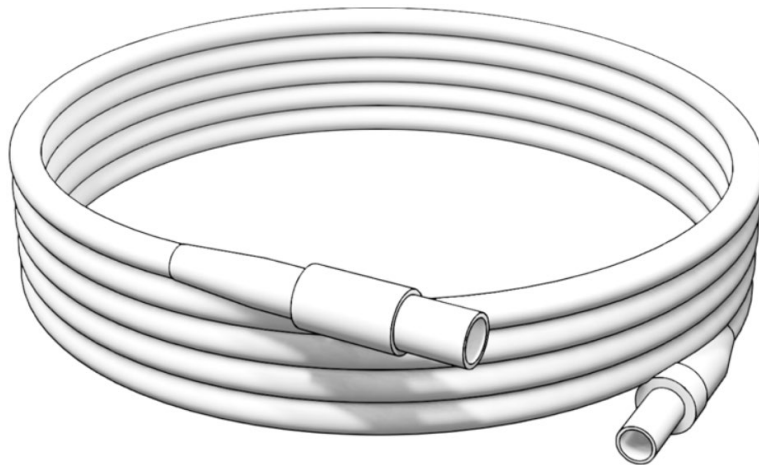English	3	Lietuvių k.	9
Česky	4	Magyar	10
Dansk	5	Nederlands	11
Deutsch	6	Norsk	12
Français	7	Português	13
Italiano	8	Suomi	14

Caution: Federal Law (USA) restricts this device to sale by or on the order of a physician.

Made in USA

DSP-0329; Rev. D
Effective date: 28 August 2025

Note: Figure is not drawn to scale.

SYMBOL LEGEND

	Authorized Representative in the European Community / Autorizovaný zástupce v Evropském společenství / Autoriseret repræsentant i EU / Autorisierte Vertretung in der Europäischen Gemeinschaft / Représentant autorisé dans la Communauté européenne / Rappresentante autorizzato per la Comunità Europea / Įgaliotasis atstovas Europos Bendrijoje / Meghatalmazott képviselő az Európai Közösségben / Gemachtigde in de Europese Gemeenschap / Autorisert representant i EU / Representante autorizado na Comunidade Europeia / Valtuutettu edustaja Euroopan yhteisössä		Catalogue Number / Katalogové číslo / Katalognummer / Katalog-Nr. / Numéro de référence / Codice catalogo / Katalogo numeris / Katalógusszám / Catalogusnummer / Katalognummer / Número de catálogo / Luettelonumero
	Caution / Upozornění / Forsigtig / Vorsicht / Attention / Attenzione / Perspėjimas / Figyelem / Let op / Forsiktig / Atenção / Huomio		Consult Operating Instructions / Prostudujte si návod k obsluze / Se betjeningsvejledningen / Bedienungsanleitung konsultieren / Consulter les consignes d'utilisation / Consultare le istruzioni per il funzionamento / Žr. naudojimo instrukciją / Ellenőrizze az üzemeltetési utasításokat / Raadpleeg de bedieningsinstructies / Les bruksanvisningen / Consultar as instruções de utilização / Katso käyttöohjeet
	Date of Manufacture / Datum výroby / Fremstillingsdato / Herstellungsdatum / Date de fabrication / Data di fabbricazione / Pagaminimo data / Gyártás dátuma / Productiedatum / Produksjonsdato / Data de fabric / Valmistuspäivämäärä		Do not use if package is damaged and consult instruction for use / Prostudujte si návod k použití. Pokud je obal poškozen, výrobek nepoužívejte / Må ikke anvendes, hvis emballagen er beskadiget, og læs brugsanvisningen / Bei beschädigter Verpackung nicht verwenden und Gebrauchsanweisung beachten / Ne pas utiliser si l'emballage est endommagé et consulter le mode d'emploi / Non utilizzare se la confezione è danneggiata e consultare le istruzioni per l'uso / Nenaudoti, jeigu pakuotė sugadinta, ir žr. naudojimo instrukcijas / Tilos felhasználni, ha a csomagolás sérült; ellenőrizze a használati utasítást / Niet gebruiken als verpakking beschadigd of geopend is en raadpleeg de gebruikersaanwijzing / Må ikke brukes dersom emballasjen er skadet, og se i bruksanvisningen / Não utilizar se a embalagem estiver danificada e consultar as instruções de utilização / Ei saa käyttää, jos pakkaus on vaurioitunut, ja katso käyttöohjeita
	Keep Away from Sunlight / Nevystavujte slunečnímu světlu / Opbevarer beskyttet mod sollys / Vor Sonnenlicht schützen / Conserver à l'abri de la lumière du soleil / Tenere lontano dalla luce solare / Laikyti atokiai nuo tiesioginės saulės šviesos / Napfénytől védve tartandó / Niet blootstellen aan zonlicht / Må beskyttes mot sollys / Manter afastado da luz solar / Säilytettävä auringonvalolta suojattuna		Keep Dry / Uchovávejte v suchu / Skal holdes tørt / Trocken lagern / Conserver au sec / Mantere asciutto / Laikyti sausoje vietoje / Szárazon tartandó / Droog bewaren / Må holdes tørr / Manter seco / Pidettävä kuivana
	Importer / Dovožce / Importør / Importeur / Hersteller / Importatore / Importuotojas / Importör / Importeur / Importør / Importador / Maahantuoja		Lot Number / Číslo šarže / Partinummer / Chargen-Nr. / Numéro de lot / Numero di lotto / Partijos numeris / Gyártásítétel-szám / Partijnummer / Partinummer / Número de lote / Eränumero
	Manufacturer / Výrobce / Producent / Hersteller / Fabricant / Produttore / Gamintojas / Gyártó / Fabrikant / Produsent / Fabricante / Valmistaja		Marking for devices entering the European market / Označení prostředků vstupujících na evropský trh / Mærkning for anordninger, der kommer ind på det europæiske marked / Markierung für Geräte, die in den europäischen Markt eintreten / Marquage des dispositifs entrant sur le marché européen / Marcatura per i dispositivi che entrano nel mercato europeo / Europos rinkoje pristatomų įrenginių žymėjimas / Az európai piacon bevezetett eszközök jelölése /

			Markering voor hulpmiddelen die op de Europese markt worden gebracht / Merking for enheter som kommer inn i det europeiske markedet / Marcação para dispositivos que entram no mercado europeu / Merkintä Euroopan markkinoilla myytävälle laitteille
	Medical Device / Zdravotnický prostředek / Medicinsk utstyr / Medizinisches Gerät / Dispositif médical / Dispositivo medico / Medicinos priemonė / Orvostechikai eszköz / Medisch hulpmiddel / Medisinsk utstyr / Dispositivo medico / Lääkinnällinen laite		Non-Sterile / Nesterilní / Ikke-steril / Nicht steril / Non stérile / Non sterile / Nesterilu / Nem steril / Niet-steriel / Usteril / Não estéril / Epästeriili
	Packaging Unit / Balení / Emballeringsenhed / Verpackungseinheit / Unité de conditionnement / Unità di confezionamento / Pakuotė / Csomagolási egység / Verpakkingseenheid / Emballasjeenhet / Unidade de embalagem / Pakkausyksikkö		Prescription Only / Pouze na lékařský předpis / Kun efter lægeordination / Verschreibungspflichtig / Uniquement sur ordonnance / Solo su prescrizione / Tik gydytojui paskyrus / Kizárolag orvosi rendelőnyre / Alleen op recept / Reseptbelagt / Apenas com receita médica / Vain lääkärin määräyksestä
	Serial Number / Sériové číslo / Serienummer / Serien-Nr. / Numéro de série / N. di serie / Serijos numeris / Sorozatszám / Serienummer / Serienummer / Número de série / Sarjanumero		Temperature Limit / Teplotní omezení / Temperaturbegrænsning / Temperaturgrenze / Limite de température / Limite di temperatura / Temperatūros apribojimas / Hőmérsékleti határérték / Temperatuurgrens / Temperaturbegrænsning / Limite de temperature / Lämpötilaraja
	Usable Length / Použitelná délka / Anvendelig længde / Verwendbare Länge / Longueur utile / Lunghezza utile / Naudojamas ilgis / Használható hossz / Bruikbare lengte / Brukslengde / Comprimento útil / Käyttöpituus		

CAUTION: Read the instructions carefully before using the product.

DEVICE DESCRIPTION

The Navistar® Cable is a component of the iCONNECT™ System. The adapter cable offers the ability to use the NAVISTAR RMT Catheter or NAVISTAR THERMOCOOL® RMT Catheter in conjunction with the iCONNECT System.

ADDITIONAL REQUIRED COMPONENTS

- iCONNECT Electronics Hub
- NAVISTAR RMT Catheter or
- NAVISTAR THERMOCOOL RMT CATHETER

INDICATIONS

The Navistar Cable is designed to connect a NAVISTAR RMT Catheter or NAVISTAR THERMOCOOL RMT CATHETER to the appropriate equipment for use with the iCONNECT System.

INTENDED PATIENT POPULATION

The intended patient population for the iCONNECT System Cable is patients undergoing diagnostic and interventional procedures in the following areas: right and left heart, and the coronary, peripheral, and neurovasculature.

INTENDED USERS

The iCONNECT System Cable should be used only by qualified medical professionals who have been thoroughly trained in its use.

CONTRAINDICATIONS

The Navistar Cable has no known contraindications.

HOW SUPPLIED

This cable is provided non-sterile. Inspect the packaging for damage prior to use. Examine the cable prior to use or reuse for any damage or defects (i.e. bent pins, kinks, nicks, crushed or elongated sections). Do not use if any damage or defects are identified.

HANDLING AND STORAGE

The product may only be used in a medical treatment facility that has been set up for the appropriate application and by trained staff. Store the product in its original packaging in a cool, dry place away from light at between 10°C and 25 °C.

AUTOClave

For steam sterilization cycles 1 through 9, place the cleaned devices in the autoclave on the bottom shelf. For 10th steam sterilization cycle: Single pouch each cable using an FDA approved pouch. Coil the cable into a loop then place it into the pouch. Do not fold the cable because doing so may break the cable. Use a pouch large enough to prevent stress to the pouch seams and to prevent excessive bending of the cable. A pouch size of 19 cm x 33 cm (7.5" x 13") or larger is recommended. Place the pouched device in a sterilizer. Do not overlap pouched devices. Sterilize the devices using the following moist heat sterilization cycle parameters:

Cycle Type	Configuration	Exposure Temperature	Pulses	Exposure Time	Minimum Drying Time
Prevacuum	Single-pouch	132°C (270°F)	4	4 minutes	20 minutes

CLEANING

The product may be contaminated after use and should be cleaned with a cloth or swab dampened in pH neutral, EPA-approved hospital grade solutions.

DISPOSAL

If damaged, dispose of the product and its residual elements or waste items in accordance with hospital regulations or local government policy.

WARNINGS

- The Navistar Cable should only be used by trained physicians.
- Dispose of product packaging in accordance with hospital, administrative and/or local government policy.
- If product damage or defects are identified, dispose of the product in accordance with hospital, administrative and/or local government policy.
- No modification of this device is allowed.
- Patient or operator injury can result from improper handling of the cable.
- Failure to abide by the above warnings might result in damage to the product, or result in serious adverse events.

CAUTIONS

- Do not immerse cable connector(s) in fluids.
- If the cable is used in the presence of electrical equipment, noise may be induced into the cable.
- Do not attempt to repair any damage. In case of doubt, discard the cable and do not use or reuse.

INSTRUCTIONS FOR USE

Refer to the iCONNECT User Guide for cable connection instructions and operation of the iCONNECT System.

TRADEMARKS

Stereotaxis, the Stereotaxis logo, iCONNECT, Genesis RMN™, and Niobe™ are trademarks or registered trademarks of Stereotaxis, Inc in the USA and other countries. All other brand names, product names, or trademarks are the property of their respective owners.

STEREOTAXIS TECHNICAL SUPPORT

For technical support, please contact Stereotaxis TeleRobotic Support Team (TST) at 1-866-269-5268 or 1-314-678-6200 or email tst@stereotaxis.com.

NOTICE TO THE USER AND/OR PATIENT

Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

DISCLAIMER OF WARRANTY AND LIMITATION OF REMEDY

STEREOTAXIS MAKES NO WARRANTIES, EXPRESS OR IMPLIED, WITH REGARD TO THE DEVICE DESCRIBED IN THIS DOCUMENT. STEREOTAXIS DISCLAIMS ALL REPRESENTATIONS OR WARRANTIES, EXPRESS OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, THE WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR USE, TITLE OR NON-INFRINGEMENT, ARISING BY STATUTE OR IN LAW, OR ARISING FROM A COURSE OF CONDUCT, COURSE OF DEALING OR USAGE OF TRADE.

STEREOTAXIS, INC. SHALL NOT BE LIABLE FOR ANY DIRECT, INCIDENTAL, OR CONSEQUENTIAL DAMAGES OTHER THAN AS EXPRESSLY PROVIDED BY SPECIFIC LAW.

STEREOTAXIS IS NOT RESPONSIBLE FOR ANY DIRECT, INCIDENTAL, OR CONSEQUENTIAL DAMAGES RESULTING FROM REUSE OF THIS DEVICE.

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©Stereotaxis, Inc. 2024
MDSS (Medical Device Safety Service GmbH)