



CS NAVIGATORTM

MAPPING CATHETER

Coronary Sinus Catheter for Deep Vein Mapping



Active Tip

2F Active Tip electrode provides bipolar signal from most distal anatomy



High-Density Signals

Tight linear electrode spacing allows for precise signal mapping and localization of arrhythmia origin




Extended Reach

Enhanced placement into the smallest of CS vasculature with a 2F active tip and 4F catheter body

The Map-iT CS Navigator, also known as the CS Navigator Active Tip Diagnostic Catheter, is a fixed curve high-density linear mapping catheter with an active tip electrode and diameter tapering to only 2F. It delivers bipolar signals from previously inaccessible areas and, with its narrow electrode spacing, enables reliable localization of the arrhythmia origin. With optimal signals, the CS Navigator ensures precise mapping of the LV summit, side branches, middle cardiac vein, and vein of Marshall. All Map-iT catheters are designed with a braid-reinforced shaft for high pushability, easy placement, and consistent handling, even during long procedures.

Specifications

French Size	4F-3.3F/2F
Length	120cm
Electrodes	10 or 20
Electrode Spacing	18-2-2-2mm, 18-2-5-2mm
Curve Type	Straight

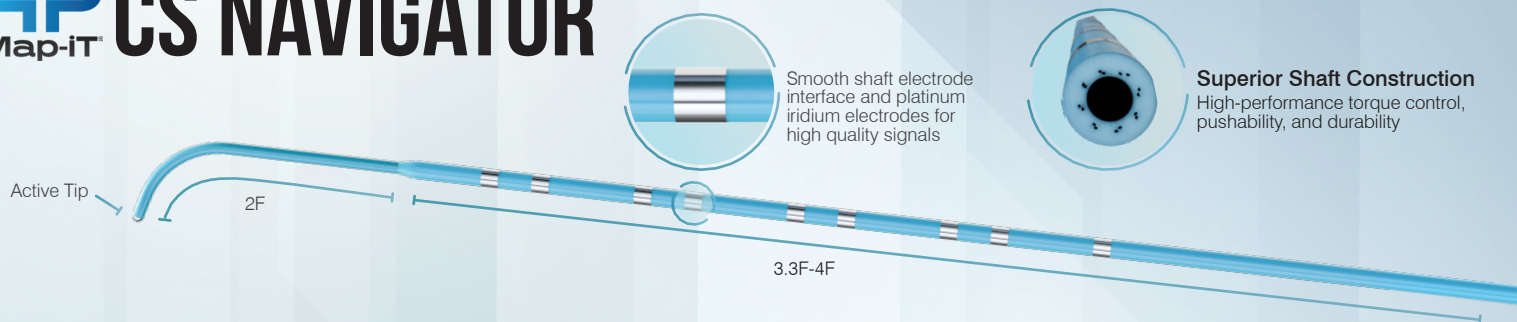


4/3.3F Fixed Diagnostic Catheter

French	4F-3.3F/2F
Ring Electrode	1mm
Tip Electrode	.66mm
Curve	Straight
Usable Length	120cm

Reference No.	Electrodes	Electrode Spacing (mm)
904222	10	18-2-2-2
904252	10	18-2-5-2
904220	20	18-2-2-2
904250	20	18-2-5-2

Connection Cable Reference No. 901957
(20 Electrode Catheter Requires 2 Cables)

Made in the USA

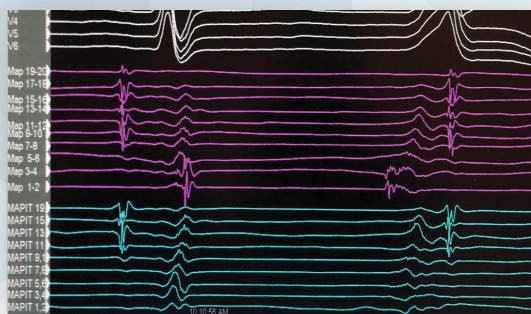


Guidewire Tip Option:

GUIDE TIP™

Coronary Sinus Catheter with .026" Guidewire Tip

Precise mapping of the LV summit, side branches, middle cardiac vein, and vein of Marshall without guide catheter support

French	4/3.3F
Ring Electrode	1mm
Tip	.026" Guidewire
Curve	Straight
Usable Length	120cm

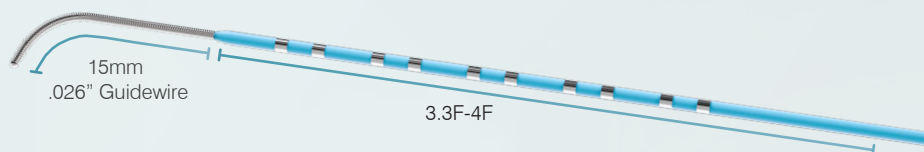

Reference No.	Electrodes	Electrode Spacing (mm)
903226	10	2-2-2
903219	10	2-5-2
903215	10	2-8-2
903217	10	5-5-5
903221	10	2-10-2
903269	20	2-2-2
903273	20	2-5-2
903256	20	2-8-2
903271	20	5-5-5
903257	20	2-10-2

Connection Cable Reference No. 901957
(20 Electrode Catheter Requires 2 Cables)

Made in the USA

To Order

✉ Map-iT@Stereotaxis.com ☎ +1 (763) 432-9004

Follow Us

✕ in EP Map-iT