



HIS RVTM

DIAGNOSTIC CATHETER

Multi-Electrode Catheter for Simultaneous HIS & RV Recording



Intentional Electrode Placement

Receive high quality electrogram data from four distal RV and six proximal HIS electrodes.



Streamlined Procedure

With a single catheter, pace the RV and record the AV node activation for a more efficient procedure.



Enhanced Access & Patient Safety

Reduce load on vasculature and preserve the opportunity for additional catheter access.


The Map-iT HIS RV, also known as the HIS-RV Diagnostic Catheter, allows physicians to use one catheter and a single access point for both HIS and RV recordings, simplifying procedures and potentially reducing the risk of HV interval prolongation. Map-iT catheters are designed with a braid-reinforced shaft for high pushability, easy placement, and consistent handling, even during long procedures.

Specifications

French Size	5F or 6F
Length	120cm
Electrodes	4 Distal, 6 Proximal
HIS to RV Electrode Spacing	20, 23, 24, 29mm
Curve Type	MPA



Smooth shaft electrode interface and platinum iridium electrodes for high quality signals

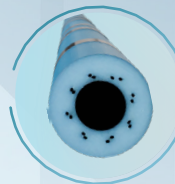

Superior Shaft Construction
High-performance torque control, pushability, and durability



Decapolar Fixed Diagnostic Catheter

Electrodes	10
Ring Electrode	1mm
Tip Electrode	2mm
Curve	MPA
Usable Length	120cm

Reference No.	French	Electrode Spacing (mm)
901579	5F	4-4-4-23-4-4-4-4-4
901580	5F	4-4-4-29-4-4-4-4-4
901583	5F	5-5-5-20-3-3-3-3-3
901584	5F	5-5-5-24-3-3-3-3-3
901823	6F	4-4-4-23-4-4-4-4-4
901824	6F	4-4-4-29-4-4-4-4-4
901827	6F	5-5-5-20-3-3-3-3-3
901828	6F	5-5-5-24-3-3-3-3-3

Connection Cable Reference No. 901957

Made in the USA



Reference Electrode Option:


HV EXPRESS[™]

Multi-Electrode Catheter for Simultaneous HIS & RV Recording with Stable Reference Electrode

Reduce load on vasculature and preserve the opportunity for additional catheter access.

Reference No.	French
950580	5
950824	6

Connection Cable Reference No. 901957

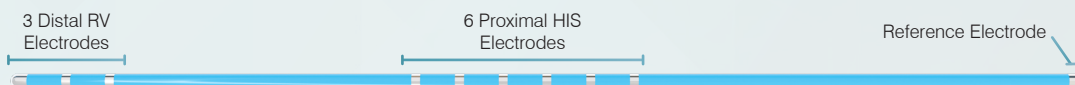

Electrodes	10
Electrode Spacing (mm)	4-4-34-4-4-4-4-4-50
Ring Electrode	1mm
Tip Electrode	2mm
Curve	MPA
Usable Length	120cm

Made in the USA

To Order

✉ Map-iT@Stereotaxis.com ☎ +1 (763) 432-9004

Follow Us

✕ in EP Map-iT